

MESA Metals: Examining Variability within Calibration Curves

Contents

Section 1: Setup	2
Section 2: Arsenic Case Study, Regression Results Replication Attempt	2
Section 3: Iterate Regression for All Arsenic	2
Section 4: Plot Arsenic Regression Lines and Intercept vs Slope	3
Section 5: All Elements Data Manipulation	4
Section 6: Plot All Elements Regression Lines and Intercept vs Slope	4

Section 1: Setup

NB a lot of the code is hidden in the PDF; investigating via visualization here.

Section 2: Arsenic Case Study, Regression Results Replication Attempt

This section reads in just one date's one calibration iteration number, for arsenic only. Then, we can fit a regression to this (the blank corrected concentrations regressed on the standard concentrations).

Printing this regression output,

```
summary(glm_fit)

##
## Call:
## glm(formula = as_o_blk_corrected ~ std_conc_ppb, family = gaussian(link = "identity"),
##      data = arsenic_As_test)
##
## Deviance Residuals:
##      1       2       3       4       5       6
## 124.99  40.84  32.65 -139.28 -177.43  118.23
##
## Coefficients:
##              Estimate Std. Error t value Pr(>|t|)
## (Intercept)   -124.99      79.70  -1.568   0.192
## std_conc_ppb  6871.96     33.92 202.606 3.56e-09 ***
## ---
## Signif. codes:  0 '***' 0.001 '**' 0.01 '*' 0.05 '.' 0.1 ' ' 1
##
## (Dispersion parameter for gaussian family taken to be 20803.51)
##
##      Null deviance: 854048044  on 5  degrees of freedom
## Residual deviance:   83214  on 4  degrees of freedom
## AIC: 80.252
##
## Number of Fisher Scoring iterations: 2
```

We can see this is not a precise replication of the first arsenic regression. This is close to the Excel sheet results but not sufficient; arsenic's 1/8/20 Calibration 1 has intercept and beta estimates of about -179.92 and 6887.72 but these show have about -124.99 and 6871.96.


I am not able to replicate the Excel result by trying different tuning approaches (is there a weight being implemented, maybe from the last sheet?). Additionally, this is not the same regression as the spectroscopy program uses at all, as per Rony's 4/1/21 email. This must be revisited and scaled across the document, but works to start looking at the patterns.

Section 3: Iterate Regression for All Arsenic


Let's iterate for each date and calibration number, for one element.

Section 4: Plot Arsenic Regression Lines and Intercept vs Slope

Then we look at these regression lines. Each unique combination of date and calibration number are their own regression line, as in the Excel.

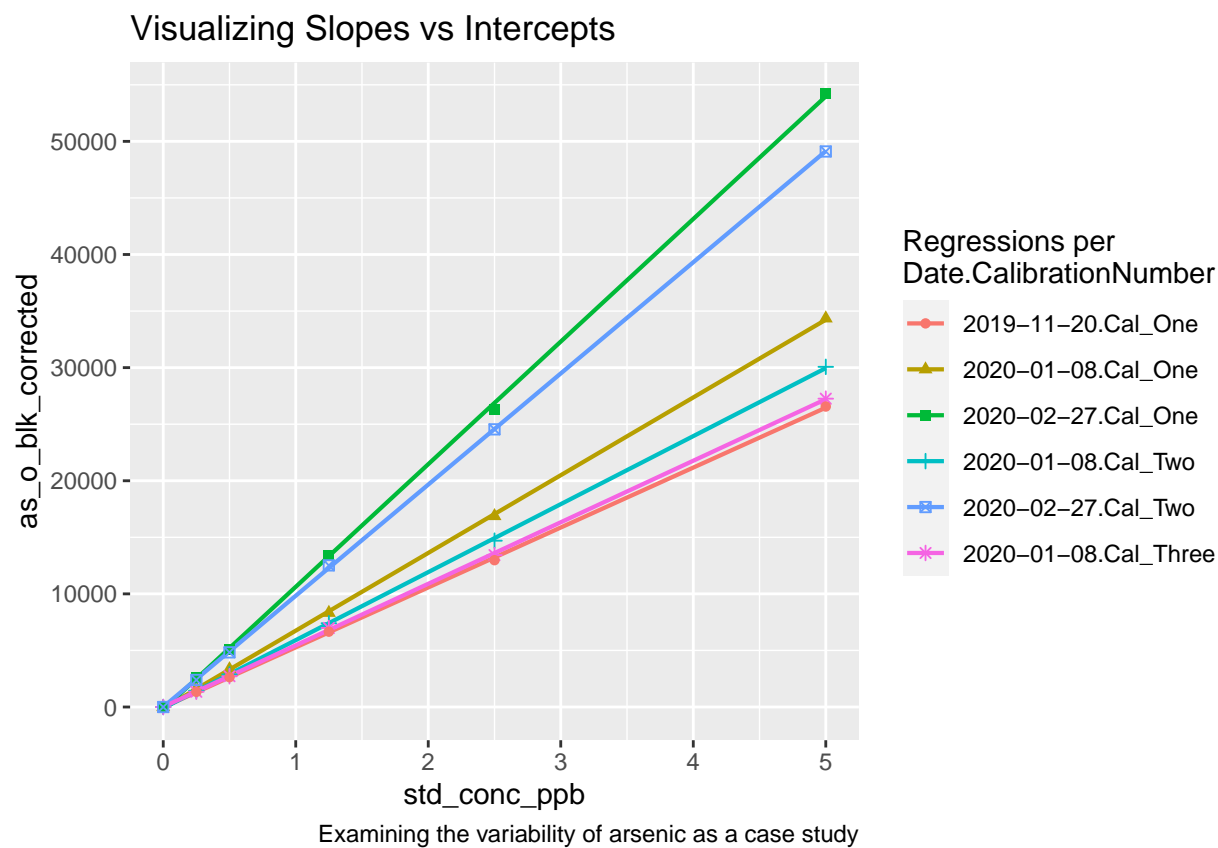

Focusing on these regression lines, let's take a look at the intercepts against the slopes.

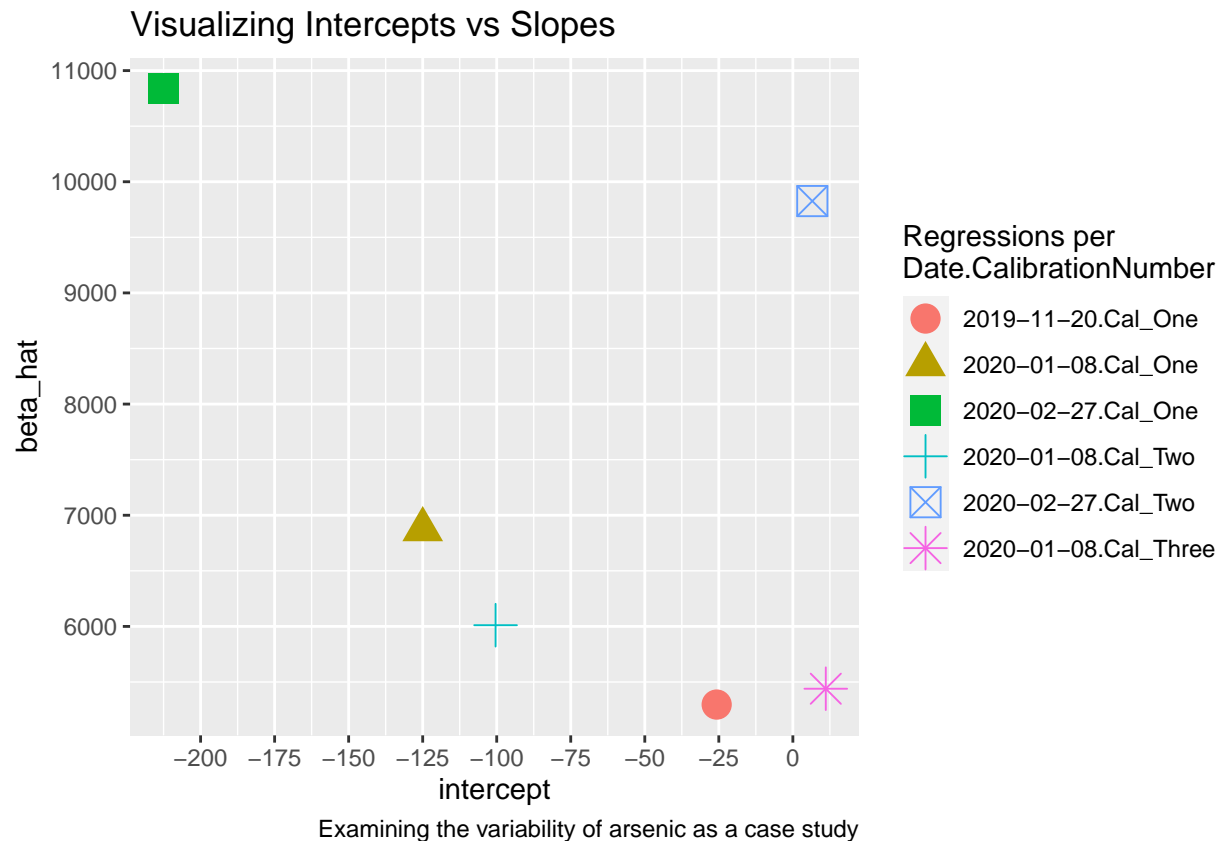

Section 5: All Elements Data Manipulation

Great. Now let's work on iterating for all elements (organized per Excel sheets within the workbook).


I read in the data.

Then, iterate the regressions for each of the (date,no_cal) groups within each element.

Section 6: Plot All Elements Regression Lines and Intercept vs Slope

Then, we're in a position to visualize the regression lines, and add the raw data to the plot just to see:


Lorem ipsum dolor sit amet, consectetur adipiscing elit. Ut purus elit, vestibulum ut, placerat ac, adipiscing vitae, felis. Curabitur dictum gravida mauris. Nam arcu libero, nonummy eget, consectetur id, vulputate a, magna. Donec vehicula augue eu neque. Pellentesque habitant morbi tristique senectus et netus et malesuada fames ac turpis egestas. Mauris ut leo. Cras viverra metus rhoncus sem. Nulla et lectus vestibulum urna fringilla ultrices. Phasellus eu tellus sit amet tortor gravida placerat. Integer sapien est, iaculis in, pretium quis, viverra ac, nunc. Praesent eget sem vel leo ultrices bibendum. Aenean faucibus. Morbi dolor nulla, malesuada eu, pulvinar at, mollis ac, nulla. Curabitur auctor semper nulla. Donec varius orci eget risus. Duis nibh mi, congue eu, accumsan eleifend, sagittis quis, diam. Duis eget orci sit amet orci dignissim rutrum.

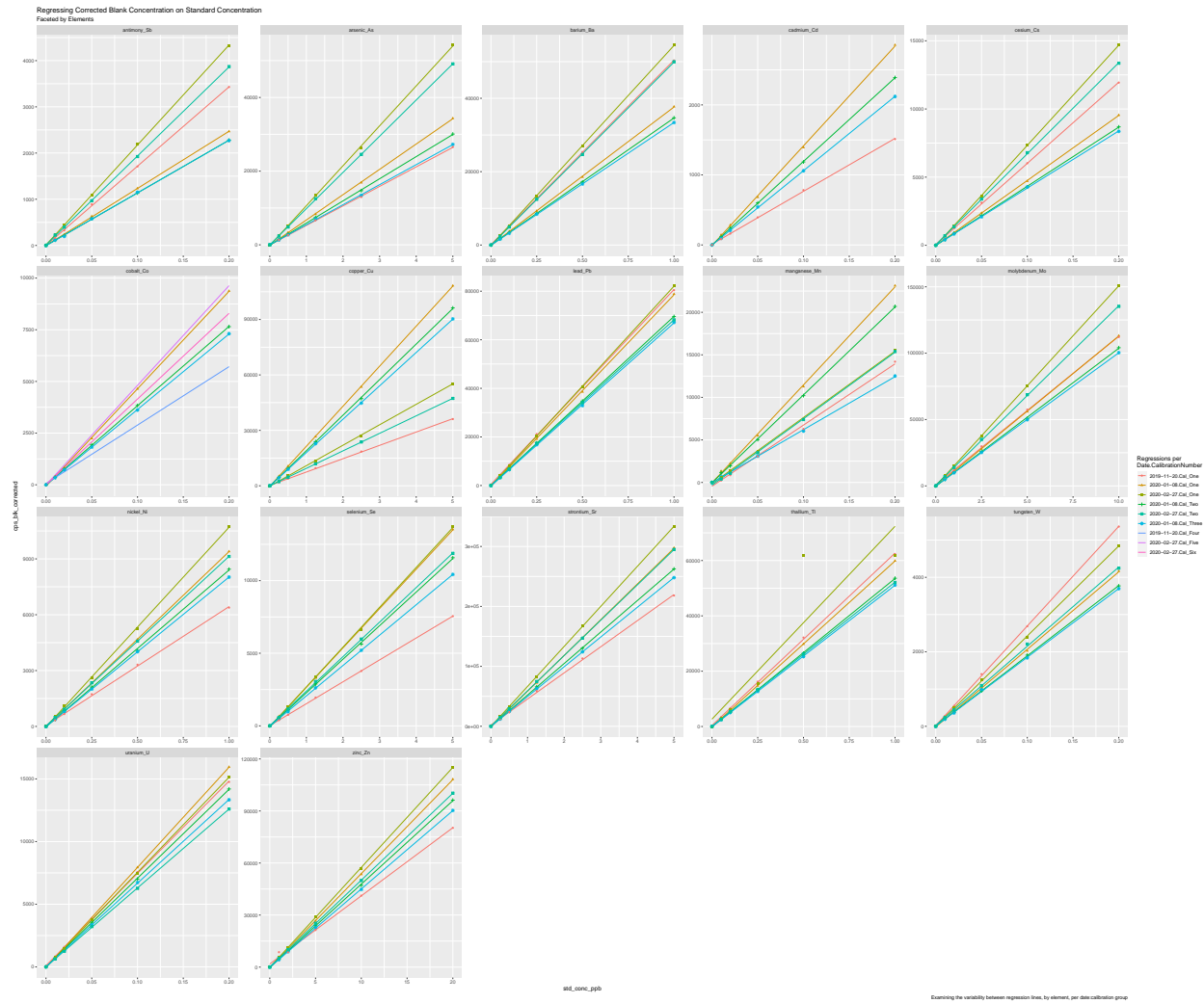

Interesting. There is some substantial variability between estimates within each element. Will think more on this. Note that thallium_Tl has a high value that is visible in the raw data (Excel too).

Like we did for arsenic, let's look at the intercepts vs slopes for these regression lines, per element also.

Visualizing intercepts vs slopes
Faceted by elements

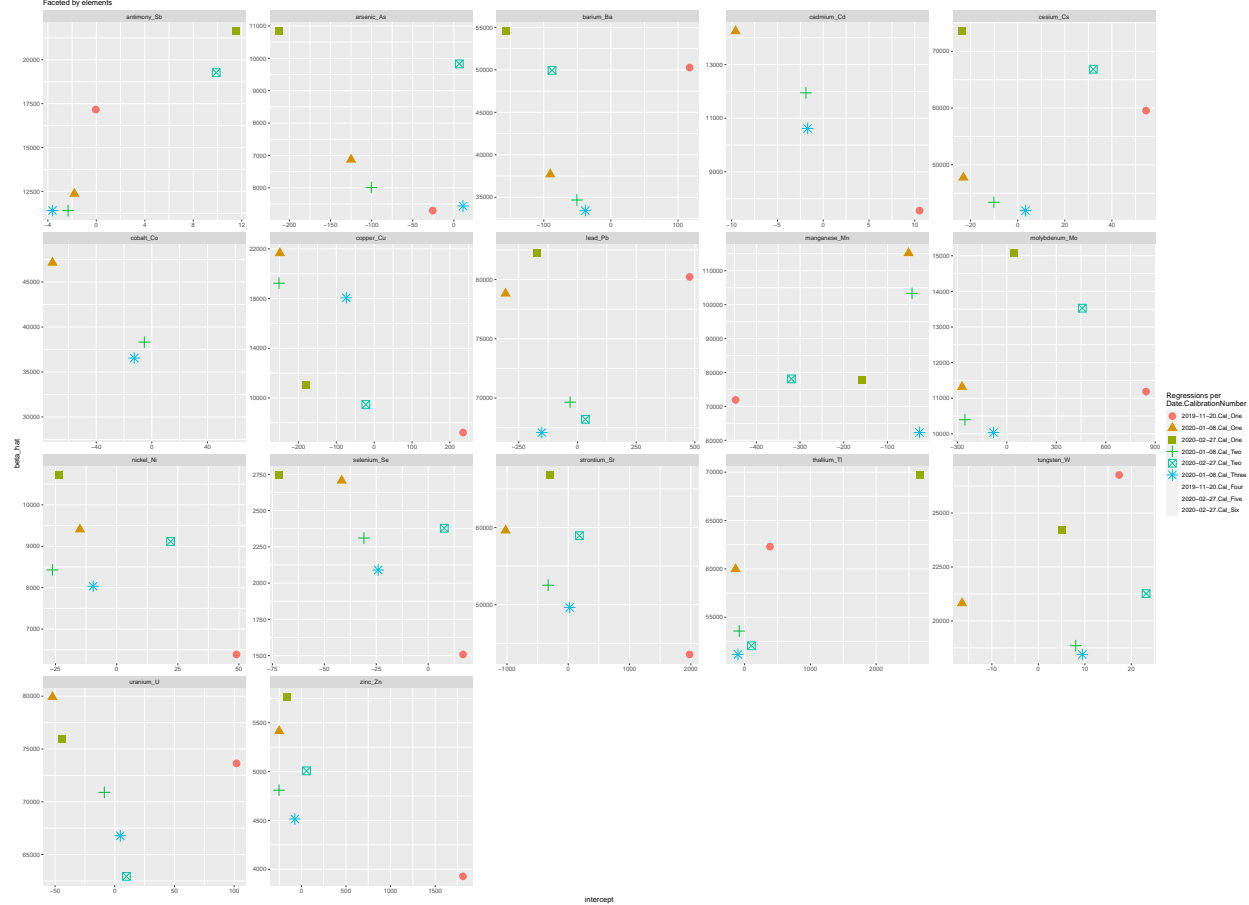

Examining the variability of intercepts vs slopes, by element, per date calibration group